





Soleo*line*, high-performance

A subtle combination of high-tech and 40 years' experience, Soleo*line* meets all the demands of modern physiotherapy.



Compact






With its large colour touch screen and SD memory card, the Soleo*line* range is extremely comfortable to use, enabling users to benefit fully from the many applications detailed in its treatment guide.

Soleo*line* is available in 3 versions: ultrasound therapy, electrostimulation and combined therapy.

Soleoline

Soleo *Sono*, innovative ultrasound therapy

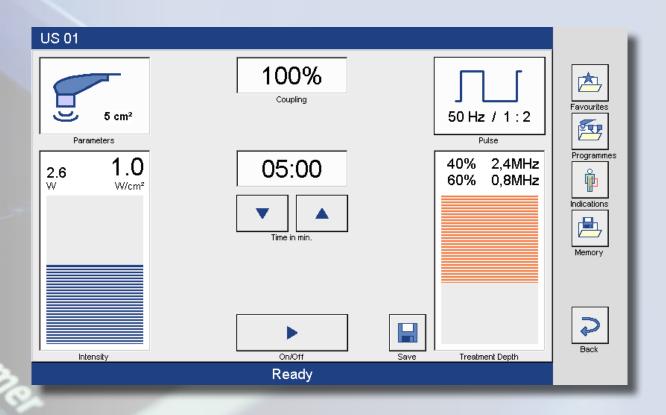



Specific programmes

The Sono comes with a very high-quality, multi-frequency (0.8 and 2.4 MHz) ergonomic head.

The preset but modifiable programmes guarantee fast and easy access to treatments.



Patent Sono Swing

With Zimmer innovation, the depth of action of the desired effect can be adjusted.

This adjustment can actually be made during the treatment for even greater precision and better therapy.



Indications

Detailed treatment recommendations are available for the main conditions affecting the locomotor system, individually stored on the SD memory card as an added bonus of the Sono.





Soleo line, user-friendly





Soleo*line* has been designed to meet all your requirements:

Speed: direct access via simple,

innovative menus

Simplicity: a large, professionally designed

touch screen


Safety: all the treatment parameters are visible

at a glance

With its intuitive software Soleo*line* is a real pleasure for practitioners to use. Its elegance and user-friendly design are features of its success

Soleoline

Soleo Galva, modern multi-channel electrostimulation

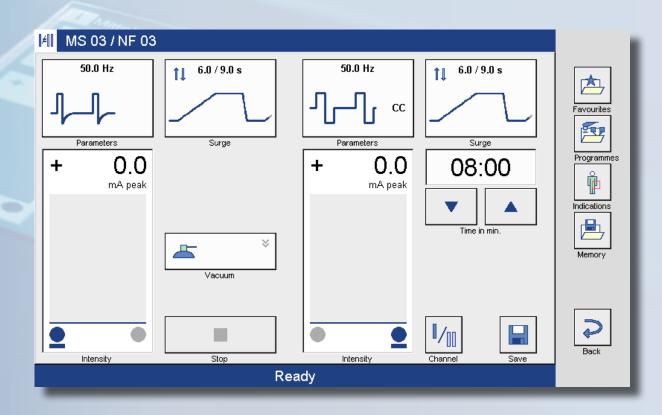




Specific programmes

"Channel" Menu

1/g Suprem

×

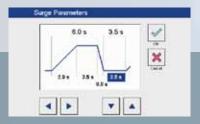

Individualised treatment





Soleo*line*, precision

Indications

Soleo line puts all the modern tools of high-performance therapy at your disposal.

40 years' experience concentrated into a Efficiency:

detailed treatment guide your modified or individualised treatments **Creativity:**

simply stored in the memory

Preference:


to your preferred programmes at all times



The exceptional performance of Soleo*line* combined with its simplicity guarantees that practitioners provide excellent treatment for the greater satisfaction of their patients. Less time spent making adjustments, more time for patients: the winning formula of the Soleo line range!

Soleolii

Soleo SonoStim, a combined device – but so simple!

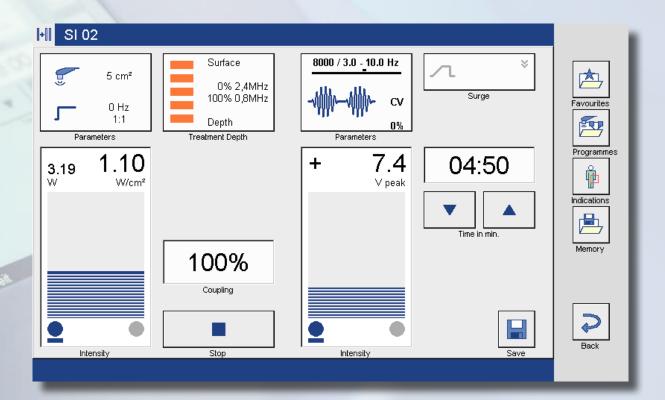



Specific programmes

The combination of Sono and Galva guarantees unparalleled flexibility and variety of treatment programmes.

In addition to ultrasound therapy and electrostimulation, Sono *Stim* offers outstanding specific combined treatments. An undeniable therapeutic bonus.



Memory

The SD memory card allows programmes and individualised 1 or 2 channel or even combined treatments to be stored.

A modern solution enabling practitioners to utilise the full capacity of the Sono *Stim* and easily monitor patients' treatment parameters.



Treatment guide

Quick selection by anatomical region gives access to more than 300 treatment recommendations for conditions affecting the locomotor system.

Myofascial pain is also catered for with the specific combined treatment programmes of the Sono *Stim*.

Soleo SonoStim







Zimmer MedizinSysteme GmbH Junkersstraße 9

Junkersstraße 9 D-89231 Neu-Ulm Tel. +49 7 31. 97 61-291 Fax +49 7 31. 97 61-299 export@zimmer.de



Specifications

Jltrasound	therapy version	

Frequencies Ultrasound head

Maximum power emission

Modes

Sono Swing patent Coupling Programmes Settings User guide

Sono

0.8 MHz: deep effect / 2.4 MHz: superficial effect 1 socket for one 5 cm² multi-frequency watertight head,

(1 cm² as an optional extra) 3 W/cm² for the 5 cm² head

Continuous and pulsed 20 Hz, 50 Hz, 100 Hz

Pulse: 1/2, 1/3, 1/5, 1/10 Selection of depth of effect Visual (2-colour) and audio 9 specific programmes

All modifiable parameters are accessible via the touch screen

More than 50 therapeutic indications

Electrostimulation version

Channels
Electrostimulation
Current types
Programmes
Specific modules
Settings

Galva

2 independent channels (programmes, intensities, parameters) Synchronised, Alternate, Separate, Sequential

Galv., LF, biphasic, MF, MF/LF, Interferential More than 100 specific programmes Muscle rehabilitation / Sports rehabilitation

All modifiable parameters are accessible via the touch screen

More than 120 therapeutic indications

Combined version

General features Programmes User guide

Sono Stim = Sono + Galva

Same features as Sono and Galva More than 110 programmes, including 5 for combined therapy

More than 300 therapeutic indications

Soleo*line*

User guide

SD memory card Accessory support Dimensions

For all versions

1 GB SD memory card / 120 memory entries + 120 favourites entries 2 supports supplied with all versions

Width 33 cm / Height 13 cm / Depth 22 cm / Weight 2 kg

Vaco equipment

Operating modes Channels Application vacuum Use Control

Dimensions

Vaco S for Galva or Sono *Stim* version Continuous and pulsed of 1 to 8 seconds

2 independent channels Adjustable from 1 to 60 kPa

Simplified by an application vacuum preset to 16 kPa $\,$

On Galva or Sono Stim screen

Width 33 cm / Height 13 cm / Depth 22 cm / Weight 3 kg

Options / Accessories

Cart
Electrodes
Remote control

1 drawer with glass tray + 1 low glass shelf

Zimmer adhesive electrodes

Wireless; with 2 channels, Intensity and Start / Stop